

t h e c o m f o r t s e r i e s ™

*My*CPAP

Presented With Compliments From  
www.MyCPAP.com



## lighten up



We're taking comfort to new heights with the Respironics Profile™ Lite Nasal Gel Mask. Softer. Lighter. More comfortable – right out of the box. The Profile Lite is also fully customizable for ultimate comfort and a personal fit. With its enhanced silent exhalation port and multiple spacer options, Profile Lite is designed to help ensure long-term compliance.

## The Profile™ Lite Mask: One-of-a-kind comfort. Two-in-one convenience.

With its softer, lighter cushion, Respironics Profile Lite Nasal Gel Mask is extremely comfortable right out of the box. The unique gel cushion and forehead pad gently rest against the patient's face for a soft, pillowy fit.

Not only is the Profile Lite a comfortable off-the-shelf mask, it's also a fully customizable mask. The simple, 10-minute customization process molds the mask to the individual's facial contours for a personalized seal. Your patients can easily customize the mask on their own, repeating the process as needed, to find the perfect fit. A personalized seal and a better fit means fewer visits and improved operating costs for you.

Most patients fit one of three sizes so there is no need to stock a large inventory of numerous sizes. If you do come across a special size requirement, other sizes can be special ordered or a standard size can be customized.

The Profile Lite consists of four parts and just two that need to be cleaned, making it easier to use and care for than other masks. Also, the Profile Lite requires less tightening to ensure a seal, thus improving patient comfort. This may lead to improved patient compliance.

### Ordering Information

	Part Number
Small Mask Only	1002850
Medium Small Mask Only	1002370
Medium Mask Only	1002371
Small Mask with Headgear	1004087
Medium Small Mask with Headgear	1004088
Medium Mask with Headgear	1004089



Special sizes available upon request.

### Replacement Parts


Deluxe Headgear	1002800
Exhalation Device	1004395
Large Gel Cushion Pad	1004291
Small Gel Cushion Pad	1004292
Extra Large Gel Cushion Pad	1008869

Covered by one or more of the following patents:  
5,540,223; DES. 377,089; 5,647,357; 5,884,624  
Other patents pending

**CAUTION:** U. S. federal law restricts this device to sale by or on the order of a physician.

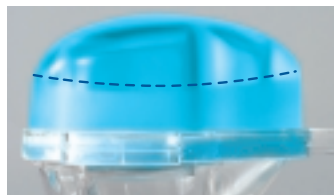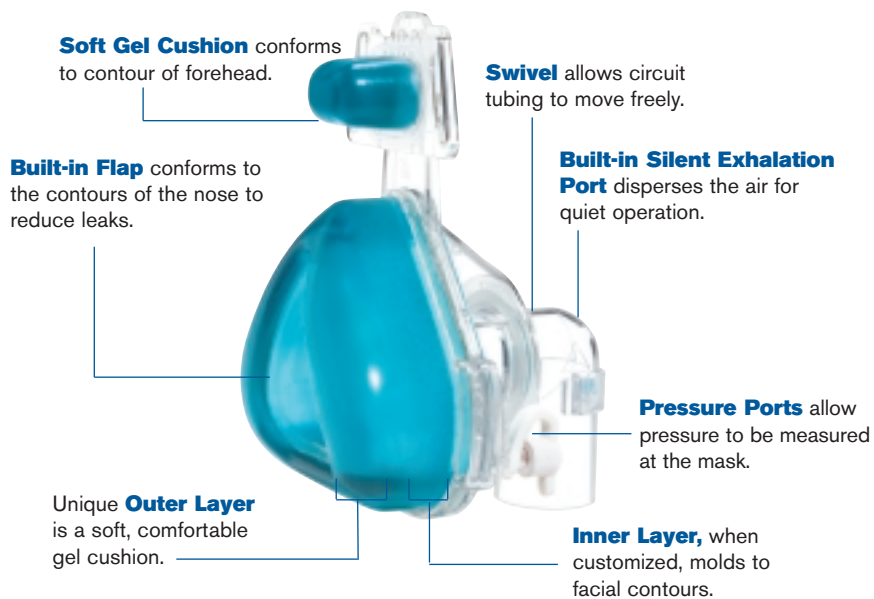

Before molding: standard contour shape.

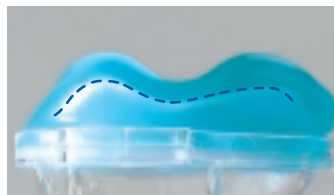

After molding: customized contour shape.

To customize Profile Lite, place the mask in water that has been boiled. Cool the mask and press it against the patient's face. The heat-absorbent cushion feels like a warm compress as it creates the custom shape.

Profile Lite comes with **Respironics** one-size-fits-most **Deluxe Headgear**. This unique headgear is adjustable in eight places allowing for a secure, comfortable fit for a wide range of patients. The Deluxe Headgear evenly distributes tension, adds stability and helps prevent mask leaks. Velcro® tabs easily adjust for a custom patient fit. The Profile Lite Nasal Mask and Deluxe Headgear disassemble easily for cleaning and replacement. The Deluxe Headgear can be used with most Respironics Masks.



Customer Service: 1-800-345-6443 • 724-387-4000  
Respironics Europe: +33-(0)1-55-60-19-80  
Respironics Asia Pacific: +852-234-342-18

Profile and Comfort Series are trademarks of Respironics, Inc.  
Velcro® is a trademark of Velcro Industries.  
Latex Free  
©2002 Respironics, Inc.

1009855 KW 03/06/02



SLEEP  
MANAGEMENT PROGRAM



CHRONIC RESPIRATORY  
MANAGEMENT PROGRAM



VENTILATION  
MANAGEMENT PROGRAM



ASTHMA, ALLERGY & SINUSITIS  
MANAGEMENT PROGRAM



INFANT  
MANAGEMENT PROGRAM