ZEST™ NASAL MASK PATIENT INSTRUCTIONS FOR USE

The Zest™ Nasal Mask is intended for single patient adult use by individuals who have been diagnosed by a physician as requiring CPAP or Bi-Level ventilator treatment.

OPERATING INSTRUCTIONS

- The operating pressure range of the mask is 3-25 cm H₂O.
- The operating temperature range of the mask is 40-104°F (5-40°C).
- Before using the mask each time:
 - Inspect it for damage. If there is any visible deterioration (cracking, tears, etc) do not use and seek replacement part(s).
 - Inspect the vent holes. Do not use if blocked.
 - Make sure the air path through the mask is clear.

Note: Failure to follow the operating instructions above may compromise the performance and safety of the mask.

WARNINGS

- This mask should be used only with CPAP or Bi-Level equipment recommended by your physician or respiratory therapist. A mask should not be used unless the CPAP or Bi-Level machine is turned on and operating properly.

Explanation of Warning: CPAP and Bi-Level machines are intended to be used with special masks which have exhalation ports to allow continuous flow of air out of the mask. When the CPAP/Bi-Level machine is turned on and functioning properly, new air from the CPAP/Bi-Level machine flushes the exhaled air out through the mask exhalation ports. However, when the CPAP/Bi-Level machine is not operating, enough fresh air will not be provided through the mask, and exhaled air may be re-breathed. Re-breathing of exhaled air for longer than several minutes can in some circumstances lead to suffocation.

- At low CPAP/Bi-Level pressures the flow through the exhalation ports may be inadequate to clear all exhaled gas from the tubing. Some re-breathing may occur.

- Do not block the gas exhaust ports. Before use, always ensure that gas is flowing out the exhaust ports.

- Do not use the mask if vomiting or feeling nauseous.

- Discontinue use if discomfort or irritation occurs, and consult your healthcare provider.

- Discontinue using the mask if there is an allergic reaction to any part of the mask. Consult your physician if this occurs.

- The Zest™ Nasal Mask is for single patient use only and cannot be reprocessed for use with different patients.

- Use with Stroke Patients: One study has shown poor compliance with CPAP therapy in stroke patients. Careful evaluation of the patient's capabilities for CPAP therapy should be exercised.

CAUTIONS

- Use the mask for its intended use only as directed in this booklet.
- US Federal Law restricts this device to sale by, or on order of, a physician.

US PATENT NUMBER 6 615 834
Other Patents Pending

FOR ADDITIONAL INFORMATION CONTACT YOUR LOCAL HOMECARE PROVIDER

FITTING YOUR MASK

- Remove plastic cover from Silicone Seal.
- Holding your mask in one hand, spread the Stretchgear™ Headgear (H) with your free hand, just enough to fit comfortably over your head.
- Lightly hold the mask against your face, covering the nose, then slowly pull the Headgear over your head.
- Attach the Headgear to the Glider™ Strap (E).
- Pull the lower straps as low as possible down the back of your neck. Gently tighten straps to remove slack. Connect to CPAP/Bi-Level machine set at prescribed pressure.
 - If necessary, tighten the lower straps to eliminate any leaks (avoid over-tightening the straps).

FITTING TIPS

- Forehead pads should rest gently on the forehead without top straps being pulled too tightly.
- If there are any leaks after gently tightening the top and lower straps, pull the mask slightly out from the face while straps are still attached. This will allow the Silicone Seal to inflate with air. Put the mask back onto the face.

CLEANING YOUR MASK

To ensure continuing safety and comfort for your Zest™ Nasal Mask:

- BEFORE FIRST USE:** Hand-wash the Stretchgear™ Headgear (H) in pure soap dissolved in lukewarm water. Do not soak for longer than 10 minutes, then rinse in fresh water and allow to dry out of direct sunlight.
- BEFORE EACH USE:** Inspect the mask for deterioration. Do not use if damaged.
- AFTER EACH USE:**
 - Disassemble the Silicone Seal (C), Foam Cushion (D) and the Headgear (H) from the Mask Base (A).
 - Wash the Silicone Seal and Mask Base in pure soap dissolved in lukewarm water. Do not soak for longer than 10 minutes.
 - Rinse thoroughly with fresh water, ensuring that all soap residue is removed.
 - Wipe the Foam Cushion with a damp cloth. (Do not immerse the Foam Cushion in water)
 - Leave all parts to dry out of direct sunlight before reassembling.

- EVERY 7 DAYS:** Hand-wash the Stretchgear™ Headgear (H) in pure soap dissolved in lukewarm water. Do not soak for longer than 10 minutes, then rinse in fresh water and allow to dry out of direct sunlight.

- Do not soak the mask for longer than 10 minutes.
- Do not wash the mask in the dishwasher.
- Do not clean the mask with products containing alcohol, anti-bacterial agents, antiseptic, bleach, chlorine or moisturizer.
- Do not store the mask in direct sunlight.

The above actions may deteriorate or damage the mask and shorten its life. If your Zest™ Nasal Mask weakens or cracks, discontinue use and seek replacement immediately.

ASSEMBLING YOUR MASK

Your Zest™ Nasal Mask is fully assembled and ready to use. However, when you need to clean your mask follow these instructions for reassembly:

- Place Silicone Seal (C) face down onto flat surface. With thick side of Foam Cushion (D) facing down, place into Silicone Seal.
- Starting with the slot at the apex of the Silicone Seal, assemble onto Mask Base (A). Push firmly together.
- Lay the Stretchgear™ Headgear (H) flat, with the loop side down. Position the assembled mask over the Headgear and attach the four straps into the corresponding slots in the Forehead Rest (F) and the Glider™ Strap (E). This can be done without undoing the Velcro™ tabs, by sliding the Headgear into the slots.

TECHNICAL SPECIFICATIONS

The Zest™ Nasal Mask and accessories are latex free. Many people find CPAP/Bi-Level therapy causes dryness of the nose and throat. Heated humidification can substantially reduce these symptoms.

- Do not soak the mask for longer than 10 minutes.
- Do not wash the mask in the dishwasher.
- Do not clean the mask with products containing alcohol, anti-bacterial agents, antiseptic, bleach, chlorine or moisturizer.
- Do not store the mask in direct sunlight.

The above actions may deteriorate or damage the mask and shorten its life. If your Zest™ Nasal Mask weakens or cracks, discontinue use and seek replacement immediately.

The Zest™ Nasal Mask and accessories are latex free. Many people find CPAP/Bi-Level therapy causes dryness of the nose and throat. Heated humidification can substantially reduce these symptoms.

OXYGEN/PRESSURE PORT CONNECTOR

If pressure readings and/or additional oxygen are required, an Oxygen/Pressure Port Connector is available (Code 900HC452).

- Connect the large end of the Connector (T) to the Swivel (G) on the mask and the small end of the Connector to the main breathing tube (R).
- Open the Port Cap(s) (P) and place the Pressure and/or Oxygen Tube (S) firmly over the Port(s) (Q).

(The two parts are identical; either part may be used for oxygen or pressure measurement. Only remove the Port Cap from parts that are being used.)


- Note:** At a fixed flow rate of supplemental oxygen flow, the inhaled oxygen concentration will vary, depending on the pressure settings, patient's breathing pattern, mask selection, and the leak rate.
- Warning:** If oxygen is used with CPAP/Bi-Level machine, the oxygen flow must be turned off when the CPAP/Bi-Level machine is not operating.

Explanation of Warning: When the CPAP/Bi-Level device is not in operation, and the oxygen flow is left on, oxygen delivered into the ventilator tubing may accumulate within the CPAP/Bi-Level machine enclosure. Oxygen accumulated in the CPAP/Bi-Level machine enclosure will create a risk of fire.

Warning: Do not smoke while using or in the vicinity of someone using oxygen with this mask.

Warning: Do not use Polyvinyl Chloride (PVC) pressure and/or oxygen tubing. Use silicone tubing instead.


Explanation of Warning: Contact with PVC can cause premature cracking of parts.

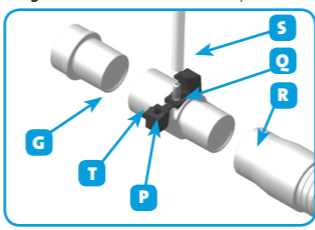EXHAUST FLOW

The Zest™ Nasal Mask has an air exhaust system that vents the expelled CO₂ through a group of fine holes in front of the mask. It is important that these holes are not blocked by any object. This controlled leak ensures exhaled CO₂ is expelled from the mask.

PRESSURE (cmH ₂ O)	3	5	7	9	11	13	15	17	19	21	23	25
FLOW (L/MIN)	17	23	28	32	36	40	44	47	49	52	55	58

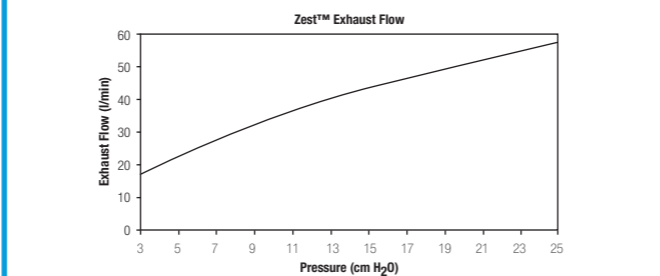RESISTANCE TO FLOW

Pressure drop through mask @ 50 l/min 0.2 cm H₂O
Pressure drop through mask @ 100 l/min 0.8 cm H₂O

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